

Prescription Misuse Prevention Kit

TimerCap®

Automatically know the last time opened to avoid medication mishaps.




TimerCap



CR Ring



Bottle

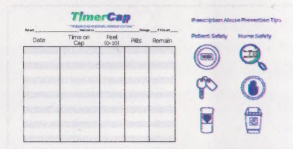

Tips and Tracker

Any questions, complaints, compliments, custom orders, insights, issues or suggestions, please contact us via e-mail at Info@TimerCap.com or visit TimerCap.com. Thank you!

TimerCap

TimerCap®

**** TAKE ONLY AS PESCIBED ****

Patient _____ Medication _____ Dosage _____ Pill Count _____

[illegible][illegible]

Prescription Misuse Prevention Tips

TimerCap®

Patient Safety



Measure

Dosage and time of medication to avoid accidental overdose.



Monitor

Your abilities before driving to avoid DUI fatalities.



Manage

Medication in containers that have label and safety information.

Home Safety



Detect

Unwanted openings by seeing the exact time it was last opened.



Deter

Misuse and theft of medication.



Dispose

Safety of unneeded medication.

TimerCap®

It's a StopWatch, Not a Clock!

No alarms and nothing to set.

HOURS MINUTES

Time Since Last Opened

Screen is blank when not on bottle and after 100 hours.

HOURS MINUTES

2:00
Time Since Last Opened

Screen shows time elapsed, **NOT** time of day.

HOURS MINUTES

SEC 1
Time Since Last Opened

Screen shows: SEC1, SEC2, then 10, 11, for the first minute.

HOURS MINUTES

99:59
Time Since Last Opened

Continues to count until 99:59, then screen will go blank.

To Start: Twist TimerCap **ON** the bottle.

To Stop: Twist TimerCap **OFF** the bottle.

If screen is blank and on the bottle, just open and close cap to start.

Really! It's that Easy!

Making TimerCap Child Restistant (Ships EZ Twist)

1.)



Put ring on top of the TimerCap.

2.)



With CR Ring on, turn upside down on a hard surface.

3.)



Click, Click
Press firmly on bottle bottom until you hear a double click.

4.)



When CR Ring is on the TimerCap properly, CR ring will not pull off.

Changing the TimerCap Battery (Battery Lasts 12 -18 Months)

1.)



Remove bottom disk and Spring.

2.)



Remove a screw holding battery clip and slide to access battery.

3.)



Replace LR-41 battery, slide metal back, & attach screw.

4.)



Snap bottom disk with Spring back in place. Done!